


Benzyl Alcohol impurities

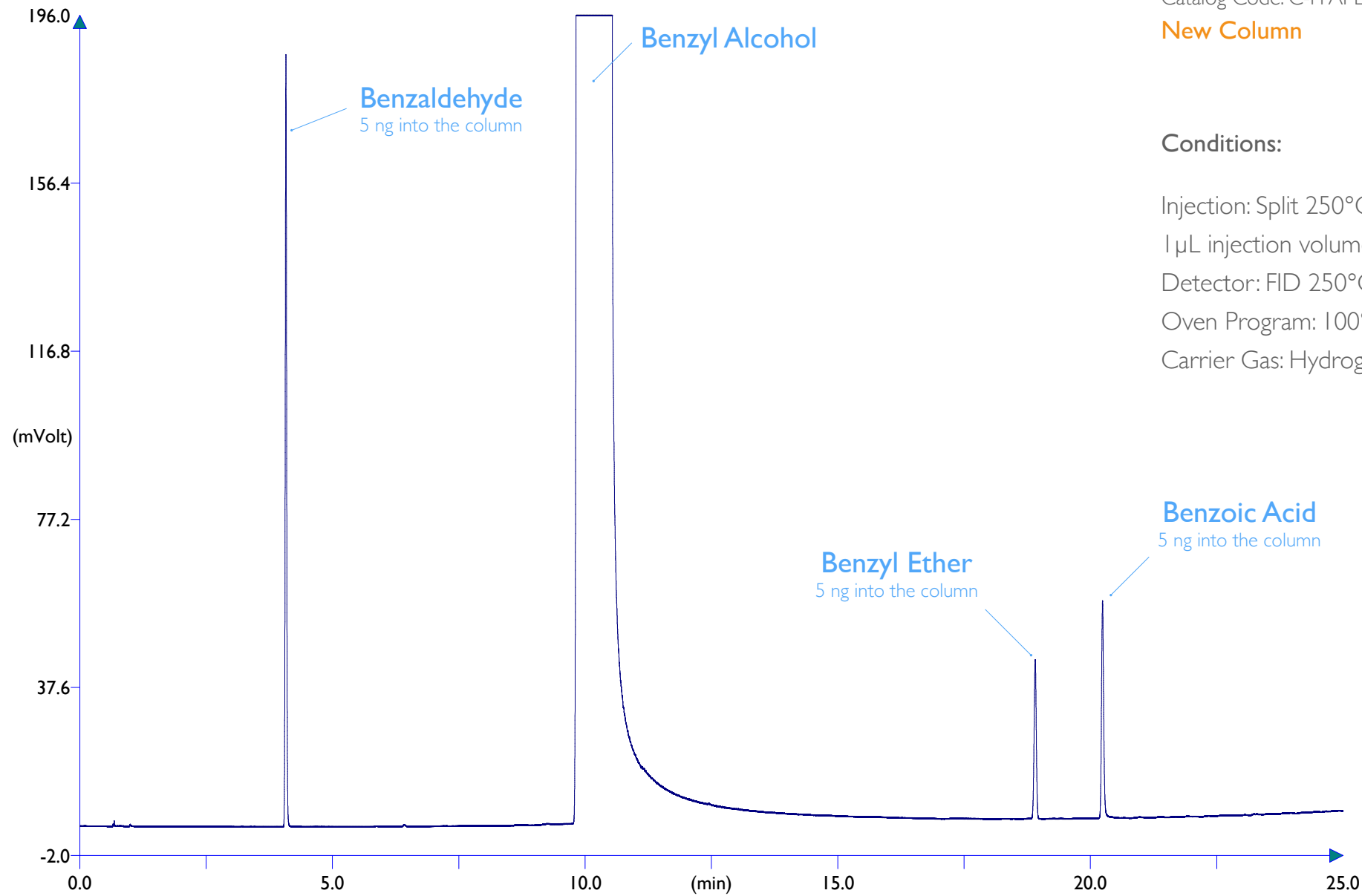

Column: **MEGA-FFAP EXT** - 0.32mm, 0.50 μ m, 30m

Catalog Code: C-FFAPEXT-032-050-30

New Column

Conditions:

Injection: Split 250°C, 80mL/min split flow,

1 μ L injection volume.

Detector: FID 250°C.

Oven Program: 100°C, 5°C/min.

Carrier Gas: Hydrogen, 80kPa.



since
1980



contact us: info@mega.mi.it

 follow us [@MEGAColumns](https://twitter.com/MEGAColumns)

www.mega.mi.it